


SAMARIT ROLLBOARDS




1. Light and convenient

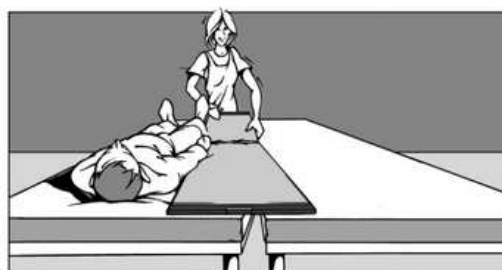

2. Position Rollbord



3. Lock brakes

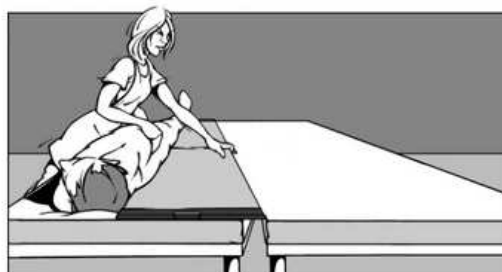

4. Turn patient to side



5. Move „Rollbord“ and ease patient back



6.1 Push patient at hip and shoulder - „One Person Push Method“



6.2 Pull patient over „Rollbord“ - „Two Person Pull Method“ USA/CA



7. Turn patient and remove „Rollbord“

Cleaning Instructions

ENGLISH

Preparation and disinfection of "SMART Rollboard"

According to recommendation and the Guidelines for Hospital Hygiene and Infection Prevention, contamination with viruses/bacteria such as MRSA, ELBL and VRE must be prevented.

We recommend the following measures:

1. Use the disposable SMART protective cover

This PE film covers the entire product and absorbs impurities, so that the Rollboard surface stays clean.

2. Wipe disinfection cont: 10-20% Isopropanol

Wipes e.g. the "Sani-Cloth " Active" or similar.

Chemical and thermal procedures are not suitable . Chlorine, concentrates and bleaches do harm materials and are used at your own risk.

FRANÇAIS

Préparation et désinfection des systèmes "Rollboards SMART"

Conformément aux directives sur l'hygiène dans le domaine hospitalier ainsi que sur la prévention des infections, toute transmission de virus ou de bactéries (SARM, ELBL, ERV, etc.) doit être évitée.

Nous recommandons les mesures suivantes :

1. Utilisation d'une housse de protection jetable SMART

Ce film de protection enveloppe l'intégralité du produit et piège les impuretés afin que la surface du Rollboard reste propre.

2. Désinfection par essuyage cont : 10-20% Isopropanol

Lingettes : p. ex. "Sani-Cloth " Active" ou similaire.

Les traitements suivants ne sont pas adaptés : procédés chimiques et thermiques, chlore, produits concentrés, agents blanchissants.

ESPAÑOL

Preparación y desinfección de los sistemas de "SMART Rollboard"

De acuerdo con las recomendaciones con las Directivas sobre Higiene hospitalaria y prevención de enfermedad des contagiosas, es preciso prevenir la transmisión de virus y bacterias, como por ejemplo, MRSA, ELBL, VRE etc.

Se recomienda adoptar las siguientes medidas:

1. Empleo de la funda protectora de un solo uso de SMART

Esta cubierta protectora envuelve todo el producto y atrapa las impurezas hasta limpiar la superficie de la Rollboard.

2. Desinfección por frotamiento cont. 10-20% Isopropanol

Con paños por ejemplo, „Sani-Cloth " Active"o similar.

Es inadecuado el uso de: procesos químicos y térmicos, cloro, concentrados o agentes blanqueadores.

