


ARTHROSCOPIC LEG HOLDER

Intended Use

Support the patients leg.

Features

- Includes one (1) Arthroscopic Leg Holder Pads (BD250-P12)
- Height Adjustment Capability
- Holds legs up to 35 inches in circumference

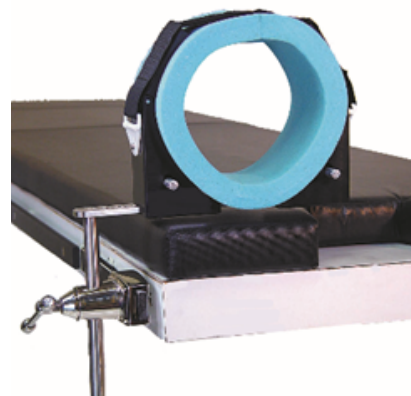

Required Equipment

- One (1) Radial Clamps such as Stainless Steel Clark Style Clamp (BD-CS)
- Arthroscopic Leg Holder Pads (BD250-P12)

⚠ Warnings

- To avoid patient injury:
 - The Leg Holder must always be used with the appropriate pad.
 - Ensure all clamps and straps are tightened before use.
 - Monitor the patient's leg for good circulation and or irritation caused by pressure during lengthy surgeries.
- The disposable Arthroscopic Leg Holder Pads are intended for one (1) time use only

Instructions

1. Loosen the strap clamp (2) that holds the pad and place a pad in the Leg Holder.
2. Insert the Leg Holder into the appropriate radial clamp on the side rail.
3. Rotate and slide the Leg Holder into position and tighten the clamp on the side rail.
4. Loosen the lateral stops (1) to allow the Leg Holder to slide into position, if necessary. Tighten the lateral stops (1) to secure the Leg Holder.
5. Place the patient's leg in the Leg Holder and wrap the pad (3) around the leg.
6. Secure the pad in place with the strap and the strap clamp (2).

Need Assistance? Please contact Technical Support (800) 804-0333 or your local representative for support